


axtair one® plus



Braden scale score

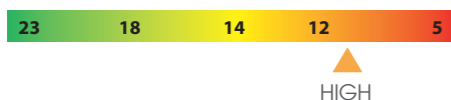

MATTRESS

- ▶ 18 removable and independent Polyurethane Ether cells of 12 cm therapeutic air
- ▶ 5 cm polyether foam base
- ▶ CPR emergency valve
- ▶ Transport plug: system for opening and sharing 2 air channels (autonomy of 8 hours)




COVERS AVAILABLE

- ▶ Upper side: stretchable, breathable, impermeable to liquids and permeable to vapour
- ▶ Lower side: non-slip base
- ▶ 2 covers available:
 - PROMUST PU
 - PROMUST CIC COVER: polyamide mesh coated polyurethane/polycarbonate with added silver ions fitted with a non-slip base



PUMP

- ▶ 1 therapeutic mode: dynamic (technology based on alternating pressure 1 cell in 2)
 - ▶ Patented system of automatic and continuous calculation of inflation pressure according to patient's morphology
 - ▶ Comfort settings: customization of comfort and seating position
 - ▶ Visual and audible alarms
 - ▶ Indication for servicing
 - ▶ Serial number on the pump's back, simplified instruction and QR code on both sides

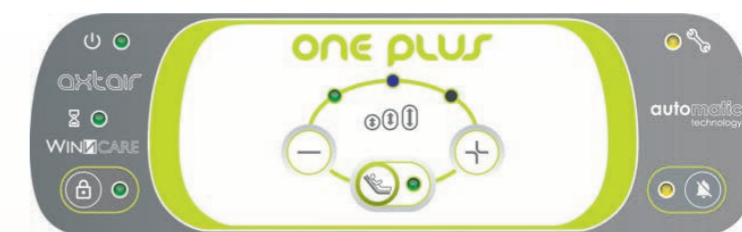

INDICATIONS

HELP IN PREVENTING PRESSURE ULCERS FOR PATIENT AT MEDIUM TO HIGH RISK OF PRESSURE ULCER DEVELOPMENT, BEDRIDDEN MORE THAN 15 HOURS A DAY.
HELP IN TREATING EXISTING PRESSURE ULCERS FROM STAGE 1 TO 2 COMBINED WITH TECHNICAL AID DEVICES.
PATIENT WEIGHT: 30 TO 150 KG/ 4.7 TO 23.6 ST

KEY BENEFITS

- ▶ **Therapeutic efficacy** for early stages of pressure ulcers
- ▶ **Easy to use** for caregivers
- ▶ **Technical autodiagnostic** by QR code
- ▶ **Ready to use** in a carrying bag
- ▶ **Environmentally-friendly**
- ▶ **Eco-designed and manufactured in France**

DESCRIPTION	REFERENCE	DIMENSIONS / WEIGHT	PATIENT WEIGHT	FIRE STANDARDS	KEY FEATURES
AXTAIR ONE® PLUS (pump + overlay mattress) Promust PU cover For 90 cm wide bed	VAXT2/ONE-P	195 x 87 x 14 cm (76.7 x 34.2 x 5.5 inches)	30 to 150 Kg	EN 597-1 and 2 GPEM D1 90 Classe D	Overlay mattress CPR emergency valve 18 removable and independent Polyurethane Ether cells of 12 cm therapeutic air (300µ) 5 cm polyether foam base Pump Dimensions: 22 x 25 x 11.5 cm Weight: 2.7 kg Technology based on alternating pressure 1 cell in 2 over a complete cycle from 9 to 14 minutes Patented system of automatic and continuous calculation of inflation pressure Comfort adjustment Visual and audible alarms. Flow > 7 l / min Electrical power supply: 220-240 V - 50 Hz Class II type BF Average consumption: 5 Wh PUMP decibels: Acoustic pressure according to NF EN3744: <35 dBA Electric and electromagnetic compliance according to the standards IEC 60601-1 Ed.3; IEC 60601-1-2
AXTAIR ONE® PLUS (pump + overlay mattress) Promust CIC cover For 90 cm wide bed	VAXT2/ONE/CIC-P	8kg			

2 years warranty

Cleaning instructions: Promust PU and CIC covers



Photos not contractually binding

FP000096/08-2018

