


# PITCHOUNE KALIN®



Bed with PITCHOUNE KALIN bed boards, white/kiwi with wooden skirting and full length side rails

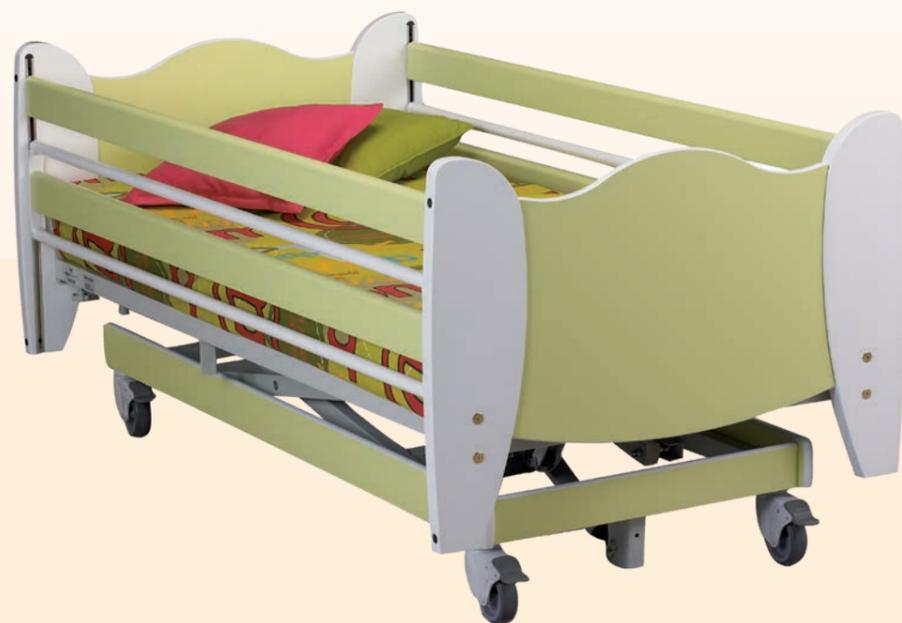

## CHILD'S BED

### FEATURES

- ▶ Non-divisible bed base 160 x 80 cm
- ▶ Base board with welded metal slat
- ▶ Equipotential plug
- ▶ 3 electric functions
- ▶ Adjustable height from 300 to 770 mm
- ▶ Electrical back rest
- ▶ Electrical leg rest with knee break by manual cranks
- ▶ On 4 x Ø 100 mm castors with independent brakes
- ▶ Centralized braking (optional)

### SAFETY

- ▶ The space between the side rails or two structure components is less than or equal to 65 mm and the space between the bed base and the bottom of the side rail is less than or equal to 65 mm
- ▶ Bed recommended for a child who is smaller than 146 cm
- ▶ Standard lockable remote control
- ▶ Optional CPR back rest

### ERGONOMICS

- ▶ Make getting up and bed-to-chair transfers easier by setting the bed base to a height of 300 mm

### STYLISH

#### 2 finishes for PITCHOUNE KALIN bed boards

- ▶ Model with white lacquered uprights with full-length wooden side rails
- ▶ Collective care model for epoxy side rails with protector for more safety



Bed with bed boards for collective care settings in ral 9016/kiwi Pitchoune Kalin, with centralized braking, epoxy side rails with protector and epoxy trapeze bar



Bed boards designed with children in mind: ral 9016/kiwi or ral 9016/Orchid Parma



PEDIA Comfort mattress in mixed viscoelastic foam, protected by a CIC cover



Option: centralized central brake system. Accessible and easy to activate

## SPECIALY DESIGNED BED

- ▶ Bed specially designed for children aged 3 to 12 years old.
- ▶ Conforms to the health & safety recommendations stipulated by AFSSaPS\* (now ANSM, France) dated 26 April 2010.

## SAFETY

- ▶ Side rails with spacing of under 65 mm.

## ERGONOMICS

- ▶ Bed base can be lowered to 300 mm.

DESCRIPTION	LPP	DIMENSIONS	BR	LEG REST	LIFTING CAPACITY	CASTORS	STANDARDS	MOTORIZATION
PITCHOUNE KALIN 3 electrical functions with high section	YES	80 x 160 cm non-divisible bed base	70°	40.3° thigh sections 14.5° adjustment with manual cranks	SWL <sup>(1)</sup> 170 kg Mattress 20 kg Patient 135 kg Accessories 15 kg	Ø 100 mm with brake without thread guard	IEC 60601-2-52 compatible with environments 1, 2, 4 and 5	LINAK IP 54
Accessories	NO	Epoxy foldaway side rails with protective cover						
	NO	PITCHOUNE KALIN for collective care						
	YES	Wooden and epoxy full-length side rails						
	YES	PITCHOUNE KALIN bed board (with spaces for full-length side rails)						
	YES	Alova PEDIA Comfort 160x80x14 cm mattress for children aged 3 to 5 and 6 to 12						
YES	HR mattress cover (included in the Pitchoune mattress pack)							

CPR back rest option: compulsory

Handset, electric and mechanical parts: 5-year warranty\*  
Boards and accessories: 1-year warranty\*  
\*Under normal conditions of use and excluding wearing parts please see instructions of use

<sup>(1)</sup> Safe working load

Photos not contractually binding

\* 26/04/2010 recommendations for users based on professional opinion and consultations with the National Health & Safety Commission for Medical Devices\*  
\*\* Depending on side rail model