



MANGAR Archimedes Bath Lift

The Archimedes Bath Lift will lower a person into a bath, on a fixed position seat, using a battery powered hand control.

EASY TO FIT

- Simple to fit – separates into three pieces, making it portable and lightweight
- Fits further back in the bath to maximise leg room and to avoid most bath handles and anti-slip coatings
- Unique 4 way quick release suckers make removing the Bath Lift simple – reducing strain on your back

SIMPLE TO MAINTAIN

- The easy to clean steel base unit has been designed with a structure that maintains strength whilst reducing weight
- No LOLER test required for two years
- Separate parts can be removed and cleaned with ease

SAFE TO USE


- Unique flexible flaps that bend on contact with an obstruction
- Bath flaps have an anti-pinch design, eliminating the risk of skin being trapped between the seat and the flap on transfer
- Only lowers a person into the bath when there is sufficient battery power to raise them back out again

IMPROVES WELL-BEING

- Bathing can reduce pain and promote a feeling of well-being
- Equipment will reduce the amount of moving and handling required from a carer
- Can be removed from the bath so enabling a more flexible environment for shared bathrooms

TECHNICAL SPECIFICATIONS

Max user weight	165 kg
Max height	45 cm
Min height	8 cm
Seat depth	44.5 cm
Total length	61 cm
Width (ex side flaps)	39.5 cm
Width of side flap	17 cm
Total product weight	10 kg
Back/headrest length	66.5 cm

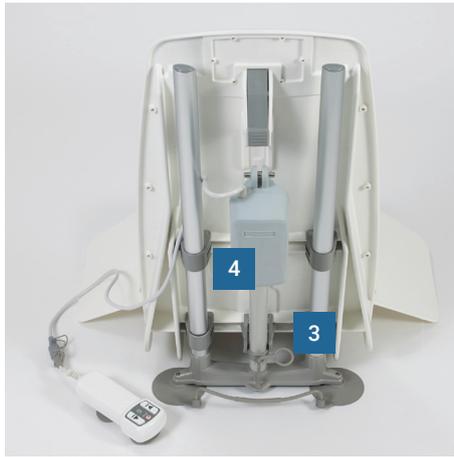

Features

1. Two anti-pinch transfer flaps
2. Sealed Hand Control with a handset which has easily visible control buttons
3. Release the transit tie by pulling up the 'O' ring and lifting the ball clear of the slot in the locating collar on the actuator.
4. Battery powered by actuator
5. Strong suckers which keeps the Archimedes in place whilst bathing
6. Safety feature which ensures there is enough charge to lift user out of bath.

Accessories

- Archimedes Spare Pack
- Archimedes Hand Control
- Archimedes Charger
- Archimedes Side Flap Protector Set



Product code: LAA3702A



For more information visit
www.mangarhealth.com
 Email: enquiries@winnicare.uk
 Phone: 0800 2800 485

